
SouthernBiotech 1



Goat Anti-Mouse IgG_{2c}, Human ads

Cat. No.	Format	Size
1079-01	Purified (UNLB)	1.0 mg
1079-02	Fluorescein (FITC)	1.0 mg
1079-04	Alkaline Phosphatase (AP)	1.0 mL
1079-05	Horseradish Peroxidase (HRP)	1.0 mL
1079-08	Biotin (BIOT)	1.0 mg
1079-09	R-phycoerythrin (PE)	0.5 mg
1079-09S	R-phycoerythrin (PE)	0.25 mg
1079-11L	Allophycocyanin (APC)	0.5 mg
1079-11S	Allophycocyanin (APC)	0.25 mg
1079-30	Alexa Fluor® 488 (AF488)	1.0 mg
1079-31	Alexa Fluor® 647 (AF647)	1.0 mg
1079-32	Alexa Fluor [®] 555 (AF555)	1.0 mg
· · · · · · · · · · · · · · · · · · ·	· · · · · · · · · · · · · · · · · · ·	·

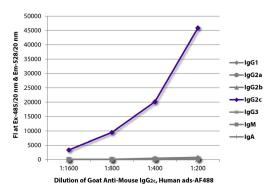

FLISA plate was coated with purified mouse IgG1, IgG2a, IgG2b, IgG_{2c}, IgG₃, IgM, and IgA. Immunoglobulins were detected with serially diluted Goat Anti-Mouse IgG2c, Human ads-AF488 (SB Cat. No. 1079-30).

Description

Specificity Reacts with the heavy chain of C57BL/6 mouse IgG_{2c}

Source Pooled antisera from goats hyperimmunized with mouse IgG_{2c}

Mouse IgG₁, IgG_{2a}, IgG_{2b}, IgG₃, IgM, and IgA; human immunoglobulins and pooled sera; may react with **Cross Adsorption**

immunoglobulins from other species

Purification Affinity chromatography on mouse IgG_{2c} covalently linked to agarose

Applications

Quality tested applications include -

ELISA 1-20 **FLISA**

Other referenced applications include -

FC ^{21,22,25} ELISPOT 5,23 IHC-PS 24 WB 25

Depletion 1,25

Working Dilutions

ELISA	AP conjugate	1:2,000 - 1:4,000
	HRP conjugate	1:4,000 - 1:8,000
	BIOT conjugate	1:5,000 – 1:20,000
FLISA	FITC and AF488 conjugates	1:200 – 1:400
	AF555 conjugate	1:100 – 1:200

PE, APC, and AF647 conjugates $\leq 1 \mu g/mL$ Since applications vary, you should determine the optimum working dilution for the product that is Other Applications

appropriate for your specific need.

For Research Use Only. Not for Diagnostic or Therapeutic Use.

Email: info@southernbiotech.com • Website: www.southernbiotech.com

Handling and Storage

- The purified (UNLB) antibody is supplied as 1.0 mg purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. No preservatives or amine-containing buffer salts added. Store at 2-8°C.
- The fluorescein (FITC) conjugate is supplied as 1.0 mg in 1.0 mL of PBS/NaN₃. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL of stock solution in 50 mM Tris/1 mM MgCl₂/50% glycerol, pH 8.0, containing NaN₃ as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL of stock solution in 50% glycerol/50% PBS, pH 7.4. No preservative added. Store at 2-8°C or long-term at -20°C.
- The biotin (BIOT) conjugate is supplied as 1.0 mg in 2.0 mL of PBS/NaN₃. Store at 2-8°C.
- The R-phycoerythrin (PE) and allophycocyanin (APC) conjugates are supplied as 0.5 mg in 1.0 mL or 0.25 mg in 0.5 mL of PBS/NaN₃ and a stabilizing agent. Store at 2-8°C. **Do not freeze!**
- The Alexa Fluor® 488 (AF488), Alexa Fluor® 555 (AF555), and Alexa Fluor® 647 (AF647) conjugates are supplied as 1.0 mg in 1.0 mL of PBS/NaN₃. Store at 2-8°C.
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

Warning

Some reagents contain sodium azide. Please refer to product specific (M)SDS.

References

- Lambert SL, Okada CY, Levy R. TCR vaccines against a murine T cell lymphoma: a primary role for antibodies of the IgG2c class in tumor protection. J Immunol. 2004;172:929-36. (ELISA. Depletion)
- 2. Žhang H, High KA, Wu Q, Yang P, Schlachterman A, Yu S, et al. Genetic analysis of the antibody response to AAV2 and factor IX. Mol Ther. 2005;11:866-74. (ELISA)
- Magalhaes JG, Fritz JH, Le Bourhis L, Sellge G, Travassos LH, Selvanantham T, et al. Nod2-dependent Th2 polarization of antigen-specific immunity. J Immunol. 2008;181:7925-35. (ELISA)
- 4. Burrell BE, Lu G, Li XC, Bishop DK. OX40 costimulation prevents allograft acceptance induced by CD40-CD40L blockade. J Immunol. 2009;182:379-90. (ELISA)
- Swanson CL, Wilson TJ, Strauch P, Colonna M, Pelanda R, Torres RM. Type I IFN enhances follicular B cell contribution to the T cell-independent antibody response. J Exp Med. 2010;207:1485-1500. (ELISA, ELISPOT)
- Kasturi SP, Skountzou I, Albrecht RA, Koutsonanos D, Hua T, Nakaya HI, et at. Programming the magnitude and persistence of antibody responses with innate immunity. Nature. 2011;470:543-7. (ELISA)
- Schmitz JM, Durham CG, Schoeb TR, Soltau TD, Wolf KJ, Tanner SM, et al. Helicobacter felis-associated gastric disease in microbiota-restricted mice. J Histochem Cytochem. 2011;59:826-41. (ELISA)
- Allan RS, Zueva E, Cammas F, Schreiber HA, Masson V, Belz GT, et al. An epigenetic silencing pathway controlling T helper 2 cell lineage commitment. Nature. 2012;487:249-53. (ELISA)
- 9. Harada Y, Tanaka S, Motomura Y, Harada Y, Ohno S, Ohno S, et al. The 3' enhancer CNS2 is a critical regulator of interleukin-4-mediated humoral immunity in follicular helper T cells. Immunity. 2012;36:188-200. (ELISA)
- Vogelzang A, McGuire HM, Liu SM, Gloss B, Mercado K, Earls P, et al. IL-21 contributes to fatal inflammatory disease in the absence of Foxp3* T regulatory cells. J Immunol. 2014;192:1404-14. (ELISA)
 Halemano K, Guo K, Heilman KJ, Barrett BS, Smith DS, Hasenkrug KJ, et al. Immunoglobulin somatic hypermutation by APOBEC3/Rfv3 during retroviral infection. Proc Natl Acad Sci
- USA. 2014;111:7759-64. (ELISA)

 12. Lam RS, O'Brien-Simpson NM, Lenzo JC, Holden JA, Brammar GC, Walsh KA, et al. Macrophage depletion abates Porphyromonas gingivalis-induced alveolar bone resorption in mice.
- 12. Lam RS, O'Brien-Simpson NM, Lenzo JC, Holden JA, Brammar GC, Walsh RA, et al. Macrophage depletion abates Porphyromonas gingivalis-induced alveolar bone resorption in mice

 J Immunol. 2014;193:2349-62. (ELISA)
- 13. Shen P, Roch T, Lampropouloù V, O'Connor RA, Stervbo U, Hilgenberg E, et al. IL-35-producing B cells are critical regulators of immunity during autoimmune and infectious diseases. Nature. 2014;507:366-70. (ELISA)
- 14. Soldemo M, Pedersen GK, Karlsson Hedestam GB. HIV-1 Env-specific memory and germinal center B cells in C57BL/6 mice. Viruses. 2014;6:3400-14. (ELISA)
- 15. Blaauboer SM, Gabrielle VD, Jin L. MPYS/STING-mediated TNF-α, not type I IFN, is essential for the mucosal adjuvant activity of (3'-5')-cyclic-di-guanosine-monophosphate in vivo. J Immunol. 2014;192:492-502. (ELISA)
- 16. Rowland SL, Riggs JM, Gilfillan S, Bugatti M, Vermi W, Kolbeck R, et al. Early, transient depletion of plasmacytoid dendritic cells ameliorates autoimmunity in a lupus model. J Exp Med. 2014;211:1977-91. (ELISA)
- Lopez-Saucedo C, Bernal-Reynaga R, Zayas-Jahuey J, Galindo-Gomez S, Shibayama M, Garcia-Galvez C, et al. CD40 ligand deficient C57BL/6 mouse is a potential surrogate model of human X-linked hyper IgM (X-HIGM) syndrome for characterizing immune responses against pathogens. Biomed Res Int. 2015;2015:679850. (ELISA)
- 18. Teo T, Her Z, Tan JJ, Lum F, Lee WW, Chan Y, et al. Caribbean and La Réunion chikungunya virus isolates differ in their capacity to induce proinflammatory Th1 and NK cell responses and acute joint pathology. J Virol. 2015;89:7955-69. (ELISA)
- Devera TŚ, Prusator DK, Joshi SK, Ballard JD, Lang ML. Immunization of mice with anthrax protective antigen limits cardiotoxicity but not hepatotoxicity following lethal toxin challenge. Toxins. 2015;7:2371-84. (ELISA)
 Blaauboer SM, Mansouri S, Tucker HR, Wang HL, Gabrielle VD, Jin L. The mucosal adjuvant cyclic di-GMP enhances antigen uptake and selectively activates pinocytosis-efficient cells
- in vivo. Elife. 2015;4:e06670. (ELISA)

 21. Bai L, Deng S, Reboulet R, Mathew R, Teyton L, Savage PB, et al. Natural killer T (NKT)-B-cell interactions promote prolonged antibody responses and long-term memory to
- pneumococcal capsular polysaccharides. Proc Natl Acad Sci USA. 2013;110:16097-102. (FC)

 22. Limon JJ, So L, Jellbauer S, Chiu H, Corado J, Sykes SM, et al. mTOR kinase inhibitors promote antibody class switching via mTORC2 inhibition. Proc Natl Acad Sci USA.
- 221. Limin 3, 30 L, Jenbauer 3, Chu H, Colaud 3, Sykes Siw, et al. INTOK kinase liminions promote animony class switching via intokez liminion. Proc Nati Acad Sci USA. 2014;111:E5076-85. (FC)

 23. Hsu H, Wu Y, Yang P, Wu Q, Job G, Chen J, et al. Overexpression of activation-induced cytidine deaminase in B cells is associated with production of highly pathogenic autoantibodies.
- J Immunol. 2007;178:5357-65. (ELISPOT)

 24. Betz BC, Jordan-Williams KL, Wang C, Kang SG, Liao J, Logan MR, et al. Batf coordinates multiple aspects of B and T cell function required for normal antibody responses. J Exp Med.
- 2010;207:933-42. (IHC-PS)
 Esser-von Bieren J, Volpe B, Kulagin M, Sutherland DB, Guiet R, Seitz A, et al. Antibody-mediated trapping of helminth larvae requires CD11b and Fcγ receptor I. J Immunol. 2015;194:1154-63. (FC, WB, Depletion)

Alexa Fluor® 488, 647, and 555 are provided under an Intellectual property license from Life Technologies Corporation. The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components for Commercial Purposes. Commercial Purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. For information on purchasing a license to this product for any other use, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com.

TB1079 13-Dec-16

Tel: 205.945.1774 • U.S. and Canada: 800.722.2255 • Fax: 205.945.8768 Email: info@southernbiotech.com • Website: www.southernbiotech.com