




Goat Anti-Mouse IgG₃, Human ads

Cat. No.	Format	Size
1100-01	Purified (UNLB)	1.0 mg
1100-02	Fluorescein (FITC)	1.0 mg
1100-03	Rhodamine (TRITC)	1.0 mg
1100-04	Alkaline Phosphatase (AP)	1.0 mL
1100-05	Horseradish Peroxidase (HRP)	1.0 mL
1100-07	Texas Red [®] (TXRD)	1.0 mg
1100-08	Biotin (BIOT)	1.0 mg
1100-09	R-phycoerythrin (PE)	0.5 mg
1100-09S	R-phycoerythrin (PE)	0.25 mg
1100-13	Spectral Red [®] (SPRD)	0.25 mg
1100-15	Cyanine 5 (CY5)	1.0 mg
1100-17	R-phycoerythrin-Cyanine 7 (PE/CY7)	0.25 mg
1100-19	Allophycocyanin-Cyanine 7 (APC/CY7)	0.25 mg
1100-30	Alexa Fluor [®] 488 (AF488)	1.0 mg
1100-31	Alexa Fluor [®] 647 (AF647)	1.0 mg
1100-32	Alexa Fluor [®] 555 (AF555)	1.0 mg

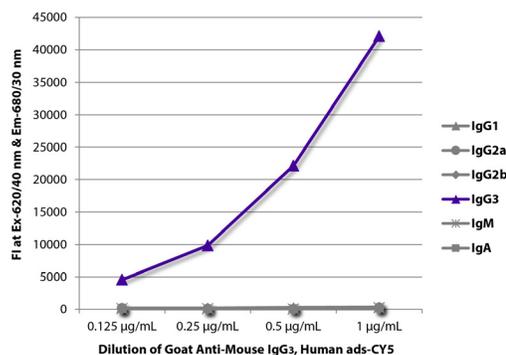

FLISA plate was coated with purified mouse IgG₁, IgG_{2a}, IgG_{2b}, IgG₃, IgM, and IgA. Immunoglobulins were detected with serially diluted Goat Anti-Mouse IgG₃, Human ads-CY5 (SB Cat. No. 1100-15).

Description

Specificity	Reacts with the heavy chain of mouse IgG ₃
Source	Pooled antisera from goats hyperimmunized with mouse IgG ₃
Cross Adsorption	Mouse IgG ₁ , IgG _{2a} , IgG _{2b} , IgM, and IgA; human immunoglobulins and pooled sera; may react with immunoglobulins from other species
Purification	Affinity chromatography on mouse IgG ₃ covalently linked to agarose

Applications

Quality tested applications include –

ELISA¹⁻⁹
 FLISA²
 FC¹⁰⁻¹⁶

Other referenced applications include –

ELISPOT^{1,3,4,10}
 IHC-FS¹⁷⁻²¹
 IHC-PS^{17,22}
 ICC²³
 WB^{1,24}

Working Dilutions

ELISA	AP conjugate	1:2,000 – 1:4,000
	HRP conjugate	1:4,000 – 1:8,000
	BIOT conjugate	1:5,000 – 1:20,000
FLISA	FITC, TRITC, TXRD, AF488, and AF555 conjugates	1:100 – 1:400
	PE, CY5, and AF647 conjugates	≤ 1 µg/mL
Flow Cytometry	FITC, BIOT, and AF488 conjugates	≤ 1 µg/10 ⁶ cells
	PE, SPRD, CY5, PE/CY7, APC/CY7, and AF647 conjugates	≤ 0.1 µg/10 ⁶ cells
	For flow cytometry, the suggested use of these reagents is in a final volume of 100 µL	
Other Applications	Since applications vary, you should determine the optimum working dilution for the product that is appropriate for your specific need.	

For Research Use Only. Not for Diagnostic or Therapeutic Use.

Handling and Storage

- The purified (UNLB) antibody is supplied as 1.0 mg purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. *No preservatives or amine-containing buffer salts added.* Store at 2-8°C.
- The fluorescein (FITC), rhodamine (TRITC), Texas Red® (TXRD), Cyanine 5 (CY5), Alexa Fluor® 488 (AF488), Alexa Fluor® 555 (AF555), and Alexa Fluor® 647 (AF647) conjugates are supplied as 1.0 mg in 1.0 mL of PBS/NaN₃. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL of stock solution in 50 mM Tris/1 mM MgCl₂/50% glycerol, pH 8.0, containing NaN₃ as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL of stock solution in 50% glycerol/50% PBS, pH 7.4. *No preservative added.* Store at 2-8°C or long-term at -20°C.
- The biotin (BIOT) conjugate is supplied as 1.0 mg in 2.0 mL of PBS/NaN₃. Store at 2-8°C.
- The R-phycoerythrin (PE) conjugate is supplied as 0.5 mg in 1.0 mL or 0.25 mg in 0.5 mL of PBS/NaN₃ and a stabilizing agent. Store at 2-8°C. **Do not freeze!**
- The Spectral Red® (SPRD) conjugate is supplied as 0.25 mg in 0.5 mL of PBS/NaN₃ and a stabilizing agent. Store at 2-8°C. **Do not freeze!**
- The R-phycoerythrin-Cyanine 7 (PE/CY7) and allophycocyanin-Cyanine 7 (APC/CY7) conjugates are supplied as 0.25 mg in 1.0 mL of PBS/NaN₃ and a stabilizing agent. Store at 2-8°C. **Do not freeze!**
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

Warning

Some reagents contain sodium azide. Please refer to product specific (M)SDS.

References

1. McLay J, Leonard E, Petersen S, Shapiro D, Greenspan NS, Schreiber JR. $\gamma 3$ gene-disrupted mice selectively deficient in the dominant IgG subclass made to bacterial polysaccharides. II. Increased susceptibility to fatal pneumococcal sepsis due to absence of anti-polysaccharide IgG3 is corrected by induction of anti-polysaccharide IgG1. *J Immunol.* 2002;168:3437-43. (ELISA, ELISPOT, WB)
2. Szmajkowski H, Smith DS, Hanson MA, Kostov Y, Lakowicz JR, Rao G. A novel method for monitoring monoclonal antibody production during cell culture. *Biotechnol Bioeng.* 2008;100:448-57. (ELISA, FLISA)
3. Hovden A, Brokstad KA, Major D, Wood J, Haaheim LR, Cox RJ. A pilot study of the immune response to whole inactivated avian influenza H7N1 virus vaccine in mice. *Influenza Other Respir Viruses.* 2009;3:21-8. (ELISA, ELISPOT)
4. Patterson HC, Kraus M, Wang D, Shahsafaei A, Henderson JM, Seagal J, et al. Cytoplasmic Iga serine/threonines fine-tune Ig α tyrosine phosphorylation and limit bone marrow plasma cell formation. *J Immunol.* 2011;187:2853-8. (ELISA, ELISPOT)
5. Förster M, Raposo B, Ekman D, Klaczkowska D, Popovic M, Nandakumar KS, et al. Genetic control of antibody production during collagen-induced arthritis development in heterogeneous stock mice. *Arthritis Rheum.* 2012;64:3594-603. (ELISA)
6. Li X, Wu J, Ptacek T, Redden DT, Brown EE, Alarcón GS, et al. Allelic-dependent expression of an activating Fc receptor on B cells enhances humoral immune responses. *Sci Transl Med.* 2013;5(216):216ra175. (ELISA)
7. Brachs S, Turqueti-Neves A, Stein M, Reimer D, Brachvogel B, Bösl M, et al. Swiprosin-1/EFhd2 limits germinal center responses and humoral type 2 immunity. *Eur J Immunol.* 2014;44:3206-19. (ELISA)
8. Pucella JN, Yen W, Kim MV, van der Veeken J, Socci ND, Naito Y, et al. miR-182 is largely dispensable for adaptive immunity: lack of correlation between expression and function. *J Immunol.* 2015;194:2635-42. (ELISA)
9. Patel PS, Kearney JF. Neonatal exposure to pneumococcal phosphorylcholine modulates the development of house dust mite allergy during adult life. *J Immunol.* 2015;194:5838-50. (ELISA)
10. Conway KL, Kuballa P, Khor B, Zhang M, Shi HN, Virgin HW, et al. ATG5 regulates plasma cell differentiation. *Autophagy.* 2013;9:528-37. (ELISPOT, FC)
11. Cook AJ, Oganessian L, Harumal P, Basten A, Brink R, Jolly CJ. Reduced switching in SCID B cells is associated with altered somatic mutation of recombined S regions. *J Immunol.* 2003;171:6556-64. (FC)
12. Martomo SA, Saribasak H, Yokoi M, Hanaoka F, Gearhart PJ. Reevaluation of the role of DNA polymerase θ in somatic hypermutation of immunoglobulin genes. *DNA Repair.* 2008;7:1603-8. (FC)
13. Johansson S, Salmon-Divon M, Johansson MH, Pickman Y, Brodin P, Kärre K, et al. Probing natural killer cell education by Ly49 receptor expression analysis and computational modelling in single MHC class I mice. *PLoS One.* 2009;4(6):e6046. (FC)
14. Brodin P, Lakshmikanth T, Johansson S, Kärre K, Höglund P. The strength of inhibitory input during education quantitatively tunes the functional responsiveness of individual natural killer cells. *Blood.* 2009;113:2434-41. (FC)
15. Pauklin S, Petersen-Mahrt SK. Progesterone inhibits activation-induced deaminase by binding to the promoter. *J Immunol.* 2009;183:1238-44. (FC)
16. Bai L, Deng S, Reboulet R, Mathew R, Teyton L, Savage PB, et al. Natural killer T (NKT)-B-cell interactions promote prolonged antibody responses and long-term memory to pneumococcal capsular polysaccharides. *Proc Natl Acad Sci USA.* 2013;110:16097-102. (FC)
17. van der Loos C. User Protocol: Practical Guide to Multiple Staining. *BioTechniques' Protocol Guide 2010.* 2009 Nov. doi: doi10.2144/000113281. (IHC-FS, IHC-PS)
18. Halstensen TS, Farstad IN, Scott H, Fausa O, Brandtzaeg P. Intraepithelial TcR α / β ⁺ lymphocytes express CD45RO more often than the TcR γ / δ ⁺ counterparts in coeliac disease. *Immunology.* 1990;71:460-6. (IHC-FS)
19. Ladinsker B, Pittschieler K. Expression of the VLA family of integrins in coeliac intestinal mucosa. *Acta Paediatr.* 1998;87:1125-30. (IHC-FS)
20. Richards HB, Satoh M, Jennette JC, Croker BP, Yoshida H, Reeves WH. Interferon- γ is required for lupus nephritis in mice treated with the hydrocarbon oil pristane. *Kidney Int.* 2001;60:2173-80. (IHC-FS)
21. Jacob N, Guo S, Mathian A, Koss MN, Gindea S, Putterman C, et al. B Cell and BAFF dependence of IFN- α -exaggerated disease in systemic lupus erythematosus-prone NZM 2328 mice. *J Immunol.* 2011;186:4984-93. (IHC-FS)
22. Chalan P, Bijzet J, Huitema MG, Kroesen B, Brouwer E, Boots AM. Expression of lectin-like transcript 1, the ligand for CD161, in rheumatoid arthritis. *PLoS One.* 2015;10(7):e0132436. (IHC-PS)
23. Haggqvist B, Hultman P. Effects of deviating the Th2-response in murine mercury-induced autoimmunity towards a Th1-response. *Clin Exp Immunol.* 2003;134:202-9. (ICC)
24. Esser-von Bieren J, Volpe B, Kulagin M, Sutherland DB, Guiet R, Seitz A, et al. Antibody-mediated trapping of helminth larvae requires CD11b and Fc γ receptor I. *J Immunol.* 2015;194:1154-63. (WB)

Texas Red® is a registered trademark of Molecular Probes, Inc.

Spectral Red® is a registered trademark of Southern Biotechnology Associates, Inc.

Spectral Red® is a PE/CY5 tandem conjugate.

Cy® is a registered trademark of GE Healthcare.

Alexa Fluor® 488, 647, and 555 are provided under an Intellectual property license from Life Technologies Corporation. The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or otherwise use this product or its components or materials made using this product or its components for Commercial Purposes. Commercial Purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. For information on purchasing a license to this product for any other use, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com.