


## Rabbit F(ab')<sub>2</sub> Anti-Sheep IgG(H+L), Human SP ads

Cat. No.	Format	Size
6016-01	Purified (UNLB)	0.5 mg
6016-02	Fluorescein (FITC)	0.5 mg
6016-04	Alkaline Phosphatase (AP)	1.0 mL
6016-05	Horseradish Peroxidase (HRP)	1.0 mL

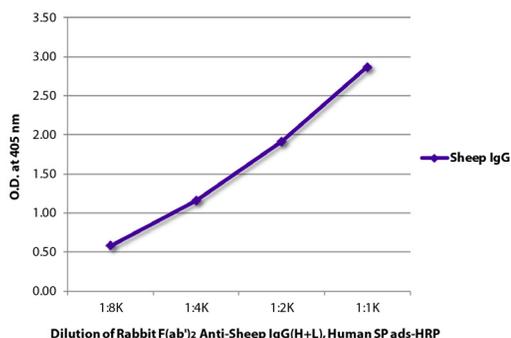

ELISA plate was coated with purified sheep IgG. Immunoglobulin was detected with Rabbit F(ab')<sub>2</sub> Anti-Sheep IgG(H+L), Human SP ads-HRP (SB Cat. No. 6016-05).

### Description

<b>Specificity</b>	Reacts with the heavy and light chains of sheep IgG
<b>Source</b>	Pepsin digest of Rabbit Anti-Sheep IgG(H+L), Human SP ads (SB Cat. No. 6156)
<b>Cross Adsorption</b>	Human serum proteins (SP); may react with immunoglobulins from other species and the light chains of other sheep immunoglobulins

### Applications

Quality tested applications include –

- ELISA <sup>1</sup>
- FLISA

### Working Dilutions

<b>ELISA</b>	AP conjugate	1:2,000 – 1:4,000
	HRP conjugate	1:4,000 – 1:8,000
<b>FLISA</b>	FITC conjugate	1:200 – 1:400
<b>Other Applications</b>	Since applications vary, you should determine the optimum working dilution for the product that is appropriate for your specific need.	

**For Research Use Only. Not for Diagnostic or Therapeutic Use.**

## Handling and Storage

---

- The purified (UNLB) antibody is supplied as 0.5 mg purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. *No preservatives or amine-containing buffer salts added.* Store at 2-8°C.
- The fluorescein (FITC) conjugate is supplied as 0.5 mg in 1.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL in a stock solution of 50 mM Tris/1 mM MgCl<sub>2</sub>/50% glycerol, pH 8.0, containing NaN<sub>3</sub> as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL in a stock solution of 50% glycerol/50% PBS, pH 7.4. No preservative added. Store at 2-8°C or long-term at -20°C.
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

## Warning

---

Some reagents contain sodium azide. Please refer to product specific SDS.

## References

---

1. Sturgill JL, Mathews J, Scherle P, Conrad DH. Glutamate signaling through the kainate receptor enhances human immunoglobulin production. *J Neuroimmunol.* 2011;233:80-9. (ELISA)