




## Goat Anti-Equine IgG(H+L)

Cat. No.	Format	Size
6040-01	Purified (UNLB)	1.0 mg
6040-02	Fluorescein (FITC)	1.0 mg
6040-04	Alkaline Phosphatase (AP)	1.0 mL
6040-05	Horseradish Peroxidase (HRP)	1.0 mL
6040-08	Biotin (BIOT)	1.0 mg
6040-09	R-phycoerythrin (PE)	0.5 mg

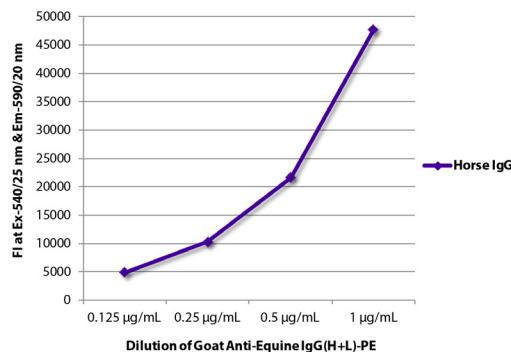

FLISA plate was coated with purified horse IgG. Immunoglobulin was detected with serially diluted Goat Anti-Equine IgG(H+L)-PE (SB Cat. No. 6040-09).

### Description

<b>Specificity</b>	Reacts with the heavy and light chains of equine IgG
<b>Source</b>	Pooled antisera from goats hyperimmunized with equine IgG
<b>Cross Adsorption</b>	None; may react with immunoglobulins from other species and the light chains of other equine immunoglobulins
<b>Purification</b>	Affinity chromatography on equine IgG covalently linked to agarose

### Applications

Quality tested applications include –

ELISA<sup>1-3</sup>  
FLISA

Other referenced applications include –

FC<sup>4,5</sup>  
ICC<sup>6</sup>  
WB<sup>7-10</sup>

### Working Dilutions

<b>ELISA</b>	AP conjugate	1:2,000 – 1:4,000
	HRP conjugate	1:4,000 – 1:8,000
	BIOT conjugate	1:5,000 – 1:20,000
<b>FLISA</b>	FITC conjugate	1:200 – 1:400
	PE conjugate	≤ 1 µg/mL
<b>Other Applications</b>	Since applications vary, you should determine the optimum working dilution for the product that is appropriate for your specific need.	

**For Research Use Only. Not for Diagnostic or Therapeutic Use.**

## Handling and Storage

---

- The purified (UNLB) antibody is supplied as 1.0 mg purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. *No preservatives or amine-containing buffer salts added.* Store at 2-8°C.
- The fluorescein (FITC) conjugate is supplied as 1.0 mg in 1.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL in a stock solution of 50 mM Tris/1 mM MgCl<sub>2</sub>/50% glycerol, pH 8.0, containing NaN<sub>3</sub> as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL in a stock solution of 50% glycerol/50% PBS, pH 7.4. *No preservative added.* Store at 2-8°C or long-term at -20°C.
- The biotin (BIOT) conjugate is supplied as 1.0 mg in 2.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- The R-phycoerythrin (PE) conjugate is supplied as 0.5 mg in 1.0 mL of PBS/NaN<sub>3</sub> and a stabilizing agent. Store at 2-8°C. **Do not freeze!**
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

## Warning

---

Some reagents contain sodium azide. Please refer to product specific SDS.

## References

---

1. Hobo S, Niwa H, Anzai T, Jones JH. Changes in serum antibody levels after vaccination for strangles and after intranasal challenge with *Streptococcus equi* subsp. *equi* in horses. *J Equine Sci.* 2010;21:33-7. (ELISA)
2. De Keyser K, Oosterlinck M, Raes E, Ducatelle R, Janssens S, Buys N. Early detection of chronic progressive lymphoedema susceptibility in Belgian draught horse stallions by means of ELISA. *Commun Agric Appl Biol Sci.* 2012;77:183-7. (ELISA)
3. De Keyser K, Berth M, Christensen N, Willaert S, Janssens S, Ducatelle R, et al. Assessment of plasma anti-elastin antibodies for use as a diagnostic aid for chronic progressive lymphoedema in Belgian Draught Horses. *Vet Immunol Immunopathol.* 2015;163:16-22. (ELISA)
4. Colson YL, Li H, Boggs SS, Patrene KD, Johnson PC, Ildstad ST. Durable mixed allogeneic chimerism and tolerance by a nonlethal radiation-based cytoreductive approach. *J Immunol.* 1996;157:2820-9. (FC)
5. Flynn J, Cox CV, Rizzo S, Foukaneli T, Rice K, Murphy M, et al. Direct binding of antithymocyte globulin to haemopoietic progenitor cells in aplastic anaemia. *Br J Haematol.* 2003;122:289-97. (FC)
6. Toscan G, Cadore GC, Pereira RC, Silva GB, Cezar AS, Sangioni LA, et al. Equine neosporosis: occurrence of antibodies against *Neospora* spp. and association between the serological status of the mares and of their offspring. *Pesq Vet Bras.* 2010;30:641-5. (ICC)
7. van 't Veer C, Mann KG. Regulation of tissue factor initiated thrombin generation by the stoichiometric inhibitors tissue factor pathway inhibitor, antithrombin-III, and heparin cofactor-II. *J Biol Chem.* 1997;272:4367-77. (WB)
8. Warner S, Hartley CA, Stevenson RA, Ficorilli N, Varrasso A, Studdert MJ, et al. Evidence that Equine rhinitis A virus VP1 is a target of neutralizing antibodies and participates directly in receptor binding. *J Virol.* 2001;75:9274-81. (WB)
9. Campbell JE, Brummel-Ziedins KE, Butenas S, Mann KG. Cellular regulation of blood coagulation: a model for venous stasis. *Blood.* 2010;116:6082-91. (WB)
10. Zhang J, Meng W, Wang C, Wu Z, Wu G, Xu Y. Expression, purification and characterization of recombinant plasminogen activator from *Gloydius brevicaudus* venom in *Escherichia coli*. *Protein Expr Purif.* 2013;91:85-90. (WB)