
# SouthernBiotech



## Rabbit Anti-Rat IgG(H+L)

Cat. No.	Format	Size
6180-01	Purified (UNLB)	1.0 mg
6180-02	Fluorescein (FITC)	1.0 mg
6180-04	Alkaline Phosphatase (AP)	1.0 mL
6180-05	Horseradish Peroxidase (HRP)	1.0 mL
6180-08	Biotin (BIOT)	1.0 mg

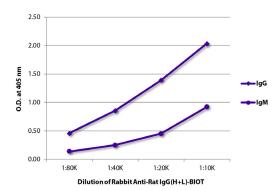

ELISA plate was coated with purified rat IgG and IgM. Immunoglobulins were detected with Rabbit Anti-Rat IgG(H+L)-BIOT (SB Cat. No. 6180-08) followed by Streptavidin-HRP (SB Cat. No. 7100-05).

## **Description**

Specificity Reacts with the heavy and light chains of rat IgG and the light chains of rat IgM

**Source** Pooled antisera from rabbits hyperimmunized with rat IgG

Cross Adsorption None; may react with immunoglobulins from other species and the light chains of other rat

immunoglobulins

**Purification** Affinity chromatography on rat IgG covalently linked to agarose

## **Applications**

Quality tested applications include -

ELISA 1

FLISA

Other referenced applications include -

IHC-FS<sup>2</sup>

IHC-PS 3-5

IP 6

WB <sup>7-10</sup>

## **Working Dilutions**

ELISA	AP conjugate	1:2,000 – 1:4,000
-------	--------------	-------------------

HRP conjugate 1:4,000 – 1:8,000 BIOT conjugate 1:5,000 – 1:20,000

**FLISA** FITC conjugate 1:200 – 1:400

Other Applications Since applications vary, you should determine the optimum working dilution for the product that is

appropriate for your specific need.

For Research Use Only. Not for Diagnostic or Therapeutic Use.

Email: info@southernbiotech.com • Website: www.southernbiotech.com

## **Handling and Storage**

- The purified (UNLB) antibody is supplied as 1.0 mg purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. No preservatives or amine-containing buffer salts added. Store at 2-8°C.
- The fluorescein (FITC) conjugate is supplied as 1.0 mg in 1.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL in a stock solution of 50 mM Tris/1 mM MgCl₂/50% glycerol, pH 8.0, containing NaN₃ as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL in a stock solution of 50% glycerol/50% PBS, pH 7.4. No preservative added. Store at 2-8°C or long-term at -20°C.
- The biotin (BIOT) conjugate is supplied as 1.0 mg in 2.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

#### Warning

Some reagents contain sodium azide. Please refer to product specific SDS.

## References

- 1. Tsuchiyama Y, Wada J, Zhang H, Morita Y, Hiragushi K, Hida K, et al. Efficacy of galectins in the amelioration of nephrotoxic serum nephritis in Wistar Kyoto rats. Kidney Int. 2000;58:1941-52. (ELISA)
- 2. Abdala-Valencia H, Earwood J, Bansal S, Jansen M, Babcock G, Garvy B, et al. Nonhematopoietic NADPH oxidase regulation of lung eosinophilia and airway hyperresponsiveness in experimentally induced asthma. Am J Physiol Lung Cell Mol Physiol. 2007;292:L1111-25. (IHC-FS)
- 3. van der Deen M, Timens W, Timmer-Bosscha H, van der Strate BW, Scheper RJ, Postma DS, et al. Reduced inflammatory response in cigarette smoke exposed Mrp1/Mdr1a/1b deficient mice. Respir Res. 2007;8:49. (IHC-PS)
- 4. Yin X, Farin HF, van Es JH, Clevers H, Langer R, Karp JM. Niche-independent high-purity cultures of Lgr5<sup>+</sup> intestinal stem cells and their progeny. Nat Methods. 2014;11:106-12. (IHC-PS)
- 5. Schaefer A, te Riet J, Ritz K, Hoogenboezem M, Anthony EC, Mul FP, et al. Actin-binding proteins differentially regulate endothelial cell stiffness, ICAM-1 function and neutrophil transmigration. J Cell Sci. 2014;127:4470-82. (IHC-PS)
- 6. Dosch DD, Ballmer-Hofer K. Transmembrane domain-mediated orientation of receptor monomers in active VEGFR-2 dimers. FASEB J. 2010;24:32-8. (IP)
- Hölscher C, Atkinson RA, Arendse B, Brown N, Myburgh E, Alber G, et al. A protective and agonistic function of IL-12p40 in mycobacterial infection. J Immunol. 2001;167:6957-66. (WB)
- 8. Gattinoni L, Klebanoff CA, Palmer DC, Wrzesinski C, Kerstann K, Yu Z, et al. Acquisition of full effector function in vitro paradoxically impairs the in vivo antitumor efficacy of adoptively transferred CD8<sup>+</sup> T cells. J Clin Invest. 2005;115:1616-26. (WB)
- 9. Shehata T, Ogawara K, Higaki K, Kimura T. Prolongation of residence time of liposome by surface-modification with mixture of hydrophilic polymers. Int J Pharm. 2008;359:272-9. (WB)
- 10. Arends F, Nowald C, Pflieger K, Boettcher K, Zahler S, Lieleg O. The biophysical properties of Basal lamina gels depend on the biochemical composition of the gel. PLoS One. 2015;10(2):e0118090. (WB)

TB6180 09-Dec-19

Email: info@southernbiotech.com • Website: www.southernbiotech.com