




Goat Anti-Mouse IgG₁, Human ads

Cat. No.	Format	Size
1070-01	Purified (UNLB)	1.0 mg
1070-02	Fluorescein (FITC)	1.0 mg
1070-03	Rhodamine (TRITC)	1.0 mg
1070-04	Alkaline Phosphatase (AP)	1.0 mL
1070-05	Horse radish Peroxidase (HRP)	1.0 mL
1070-07	Texas Red [®] (TXRD)	1.0 mg
1070-08	Biotin (BIOT)	1.0 mg
1070-09	R-phycoerythrin (PE)	0.5 mg
1070-09S	R-phycoerythrin (PE)	0.25 mg
1070-11L	Allophycocyanin (APC)	0.5 mg
1070-11S	Allophycocyanin (APC)	0.25 mg
1070-15	Cyanine 5 (CY5)	1.0 mg
1070-19	Allophycocyanin-Cyanine 7 (APC/CY7)	0.25 mg
1070-30	Alexa Fluor [®] 488 (AF488)	1.0 mg
1070-31	Alexa Fluor [®] 647 (AF647)	1.0 mg
1070-32	Alexa Fluor [®] 555 (AF555)	1.0 mg

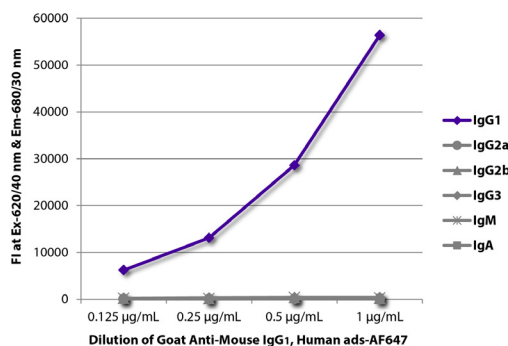

FLISA plate was coated with purified mouse IgG₁, IgG_{2a}, IgG_{2b}, IgG₃, IgM, and IgA. Immunoglobulins were detected with serially diluted Goat Anti-Mouse IgG₁, Human ads-AF647 (SB Cat. No. 1070-31).

Description

Specificity	Reacts with the heavy chain of mouse IgG ₁
Source	Pooled antisera from goats hyperimmunized with mouse IgG ₁
Cross Adsorption	Mouse IgG _{2a} , IgG _{2b} , IgG ₃ , IgM, and IgA; human immunoglobulins and pooled sera; may react with immunoglobulins from other species
Purification	Affinity chromatography on mouse IgG ₁ covalently linked to agarose

Applications

Quality tested applications include –

ELISA¹⁻⁸
 FLISA
 FC^{6,9-15}

Other referenced applications include –

ELISPOT^{4,6,7}
 IHC-FS^{3,9,16,17}
 IHC-PS^{1,16,18,19}
 IHC-WM²⁵
 ICC^{14,20,21}
 EM¹⁷
 WB^{2,5,13,22,23}
 Multiplex²⁶
 SPR²⁴

Working Dilutions

ELISA	AP conjugate	1:2,000 – 1:4,000
	HRP conjugate	1:4,000 – 1:8,000
	BIOT conjugate	1:5,000 – 1:20,000

FLISA	FITC, TRITC, TXRD, AF488, and AF555 conjugates	1:100 – 1:400
	PE, APC, CY5, and AF647 conjugates	≤ 1 µg/mL

Flow Cytometry	FITC, BIOT, and AF488 conjugates	≤ 1 µg/10 ⁶ cells
	PE, APC, CY5, APC/CY7, and AF647 conjugates	≤ 0.1 µg/10 ⁶ cells

For flow cytometry, the suggested use of these reagents is in a final volume of 100 µL

Other Applications Since applications vary, you should determine the optimum working dilution for the product that is appropriate for your specific need.

For Research Use Only. Not for Diagnostic or Therapeutic Use.

Handling and Storage

- The purified (UNLB) antibody is supplied as 1.0 mg purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. *No preservatives or amine-containing buffer salts added.* Store at 2-8°C.
- The fluorescein (FITC), rhodamine (TRITC), Texas Red® (TXRD), Cyanine 5 (CY5), Alexa Fluor® 488 (AF488), Alexa Fluor® 555 (AF555), and Alexa Fluor® 647 (AF647) conjugates are supplied as 1.0 mg in 1.0 mL of PBS/NaN₃. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL in a stock solution of 50 mM Tris/1 mM MgCl₂/50% glycerol, pH 8.0, containing NaN₃ as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL in a stock solution of 50% glycerol/50% PBS, pH 7.4. No preservative added. Store at 2-8°C or long-term at -20°C.
- The biotin (BIOT) conjugate is supplied as 1.0 mg in 2.0 mL of PBS/NaN₃. Store at 2-8°C.
- The R-phycoerythrin (PE) and allophycocyanin (APC) conjugates are supplied as 0.5 mg in 1.0 mL or 0.25 mg in 0.5 mL of PBS/NaN₃ and a stabilizing agent. Store at 2-8°C. **Do not freeze!**
- The allophycocyanin-Cyanine 7 (APC/CY7) conjugate is supplied as 0.25 mg in 1.0 mL of PBS/NaN₃ and a stabilizing agent. Store at 2-8°C. **Do not freeze!**
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

Warning

Some reagents contain sodium azide. Please refer to product specific SDS.

References

1. Puttner C, Diamond B. Immunization with a peptide surrogate for double-stranded DNA (dsDNA) induces autoantibody production and renal immunoglobulin deposition. *J Exp Med.* 1998;188:29-38. (ELISA, IHC-PS)
2. Martin RM, Brady JL, Lew AM. The need for IgG2c specific antiserum when isotyping antibodies from C57BL/6 and NOD mice. *J Immunol Methods.* 1998;212:187-92. (ELISA, WB)
3. Richards HB, Satoh M, Jennette JC, Croker BP, Yoshida H, Reeves WH. Interferon- γ is required for lupus nephritis in mice treated with the hydrocarbon oil pristane. *Kidney Int.* 2001;60:2173-80. (ELISA, IHC-FS)
4. Sanjo H, Hikida M, Aiba Y, Mori Y, Hatano N, Ogata M, et al. Extracellular signal-regulated protein kinase 2 is required for efficient generation of B cells bearing antigen-specific immunoglobulin G. *Mol Cell Biol.* 2007;27:1236-46. (ELISA, ELISPOT)
5. Yamazaki T, Nagashima M, Ninomiya D, Arai Y, Teshima Y, Fujimoto A, et al. Passive immune-prophylaxis against influenza virus infection by the expression of neutralizing anti-hemagglutinin monoclonal antibodies from plasmids. *Jpn J Infect Dis.* 2011;64:40-9. (ELISA, WB)
6. Conway KL, Kuballa P, Khor B, Zhang M, Shi HN, Virgin HW, et al. ATG5 regulates plasma cell differentiation. *Autophagy.* 2013;9:528-37. (ELISA, ELISPOT, FC)
7. Koutsouanos DG, Esser ES, McMaster SR, Kalluri P, Lee J, Prausnitz MR, et al. Enhanced immune responses by skin vaccination with influenza subunit vaccine in young hosts. *Vaccine.* 2015;33:4675-82. (ELISA, ELISPOT)
8. Dabaghian M, Latifi AM, Tebianian M, Dabaghian F, Ebrahimi SM. A truncated C-terminal fragment of Mycobacterium tuberculosis HSP70 enhances cell-mediated immune response and longevity of the total IgG to influenza A virus M2e protein in mice. *Antiviral Res.* 2015;120:23-31. (ELISA)
9. Farstad IN, Carlsen H, Morton HC, Brandtzaeg P. Immunoglobulin A cell distribution in the human small intestine: phenotypic and functional characteristics. *Immunology.* 2000;101:354-63. (FC, IHC-FS)
10. Zurney J, Howard KE, Sherry B. Basal expression levels of IFNAR and Jak-STAT components are determinants of cell-type-specific differences in cardiac antiviral responses. *J Virol.* 2007;81:13668-80. (FC)
11. Pauklin S, Petersen-Mahrt SK. Progesterone inhibits activation-induced deaminase by binding to the promoter. *J Immunol.* 2009;183:1238-44. (FC)
12. Bai L, Deng S, Reboulet R, Mathew R, Teyton L, Savage PB, et al. Natural killer T (NKT)-B-cell interactions promote prolonged antibody responses and long-term memory to pneumococcal capsular polysaccharides. *Proc Natl Acad Sci USA.* 2013;110:16097-102. (FC)
13. Sutton KM, Hu T, Wu Z, Siklodi B, Vervelde L, Kaiser P. The functions of the avian receptor activator of NF- κ B ligand (RANKL) and its receptors, RANK and osteoprotegerin, are evolutionarily conserved. *Dev Comp Immunol.* 2015;51:170-84. (FC, WB)
14. Teocchi MA, Ramalho VD, Abramczuk BM, D'Souza-Li L, Santos Vilela MM. BTK mutations selectively regulate BTK expression and upregulate monocyte XBP1 mRNA in XLA patients. *Immun Inflamm Dis.* 2015;3:171-81. (FC, ICC)
15. Imamura M, Shook D, Kamiya T, Shimasaki N, Chai SM, Coustan-Smith E, et al. Autonomous growth and increased cytotoxicity of natural killer cells expressing membrane-bound interleukin-15. *Blood.* 2014;124:1081-8. (FC)
16. van der Loos C. User Protocol: Practical Guide to Multiple Staining. *BioTechniques' Protocol Guide 2010.* 2009 Nov. doi: doi10.2144/000113281. (IHC-FS, IHC-PS)
17. Raats CJ, van den Born J, Bakker MA, Oppers-Walgreen B, Pisa BJ, Dijkman HB, et al. Expression of agrin, dystroglycan, and utrophin in normal renal tissue and in experimental glomerulopathies. *Am J Pathol.* 2000;156:1749-65. (IHC-FS, EM)
18. Huizinga R, van der Star BJ, Kipp M, Jong R, Gerritsen W, Clamer T, et al. Phagocytosis of neuronal debris by microglia is associated with neuronal damage in multiple sclerosis. *Glia.* 2012;60:422-31. (IHC-PS)
19. Neufeld T, Ludwig B, Barkai U, Weir GC, Colton CK, Evron Y, et al. The efficacy of an immunisolating membrane system for islet xenotransplantation in minipigs. *PLoS One.* 2013;8(8):e70150. (IHC-PS)
20. King MC, Raposo G, Lemmon MA. Inhibition of nuclear import and cell-cycle progression by mutated forms of the dynamin-like GTPase MxB. *Proc Natl Acad Sci USA.* 2004;101:8957-62. (ICC)
21. Schindler YL, Garske KM, Wang J, Firulli BA, Firulli AB, Poss KD, et al. Hand2 elevates cardiomyocyte production during zebrafish heart development and regeneration. *Development.* 2014;141:3112-22. (ICC)
22. Konitsiotis AD, Jovanović B, Ciepla P, Spitaler M, Lanyon-Hogg T, Tate EW, et al. Topological analysis of Hedgehog acyltransferase, a multipalmitoylated transmembrane protein. *J Biol Chem.* 2015;290:3293-307. (WB)
23. Esser-von Bieren J, Volpe B, Kulagin M, Sutherland DB, Guet R, Seitz A, et al. Antibody-mediated trapping of helminth larvae requires CD11b and Fc γ receptor I. *J Immunol.* 2015;194:1154-63. (WB)
24. Sun Y, Peng I, Senger K, Hamidzadeh K, Reichelt M, Baca M, et al. Critical role of activation induced cytidine deaminase in experimental autoimmune encephalomyelitis. *Autoimmunity.* 2013;46:157-67. (Surface Plasmon Resonance)
25. Zhu D, Fang Y, Gao K, Shen J, Zhong TP, Li F. Vegfa impacts early myocardium development in zebrafish. *Int J Mol Sci.* 2017;18:444. (IHC-WM)
26. Brummelman J, Helm K, Hamstra H, van der Ley P, Boog CJ, Han WG, et al. Modulation of the CD4⁺ T cell response after acellular pertussis vaccination in the presence of TLR4 ligation. *Vaccine.* 2015;33:1483-91. (Multiplex)

Texas Red® is a registered trademark of Molecular Probes, Inc.

Cy® is a registered trademark of GE Healthcare.

Alexa Fluor® 488, 647, and 555 are provided under an Intellectual property license from Life Technologies Corporation. The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or otherwise use this product or its components or materials made using this product or its components for Commercial Purposes. Commercial Purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. For information on purchasing a license to this product for any other use, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com.