




## Goat Anti-Bovine IgG(H+L)

Cat. No.	Format	Size
6030-01	Purified (UNLB)	1.0 mg
6030-02	Fluorescein (FITC)	1.0 mg
6030-04	Alkaline Phosphatase (AP)	1.0 mL
6030-05	Horseradish Peroxidase (HRP)	1.0 mL
6030-08	Biotin (BIOT)	1.0 mg

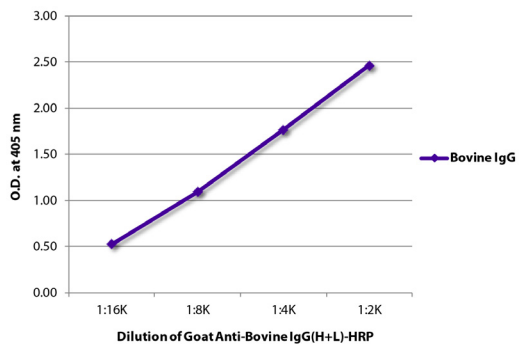

ELISA plate was coated with purified bovine IgG. Immunoglobulin was detected with Goat Anti-Bovine IgG(H+L)-HRP (SB Cat. No. 6030-05).

### Description

<b>Specificity</b>	Reacts with the heavy and light chains of bovine IgG
<b>Source</b>	Pooled antisera from goats hyperimmunized with bovine IgG
<b>Cross Adsorption</b>	None; may react with immunoglobulins from other species and the light chains of other bovine immunoglobulins
<b>Purification</b>	Affinity chromatography on bovine IgG covalently linked to agarose

### Applications

Quality tested applications include –

ELISA <sup>1-6</sup>  
FLISA

Other referenced applications include –

FC <sup>7</sup>  
ICC <sup>6,8</sup>  
IP <sup>9</sup>  
WB <sup>2,10</sup>  
Purification <sup>1</sup>  
SPR <sup>11</sup>

### Working Dilutions

<b>ELISA</b>	AP conjugate	1:2,000 – 1:4,000
	HRP conjugate	1:4,000 – 1:8,000
	BIOT conjugate	1:5,000 – 1:20,000
<b>FLISA</b>	FITC conjugate	1:200 – 1:400
<b>Other Applications</b>	Since applications vary, you should determine the optimum working dilution for the product that is appropriate for your specific need.	

**For Research Use Only. Not for Diagnostic or Therapeutic Use.**

**Corporate Offices:** 160 Oxmoor Blvd • Birmingham, AL 35209 • USA **Mailing Address:** P.O. Box 26221 • Birmingham, AL 35260 • USA

**Tel:** 205.945.1774 • U.S. and Canada: 800.722.2255 • **Fax:** 205.945.8768

**Email:** [info@southernbiotech.com](mailto:info@southernbiotech.com) • **Website:** [www.southernbiotech.com](http://www.southernbiotech.com)

## Handling and Storage

---

- The purified (UNLB) antibody is supplied as 1.0 mg purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. *No preservatives or amine-containing buffer salts added.* Store at 2-8°C.
- The fluorescein (FITC) conjugate is supplied as 1.0 mg in 1.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL in a stock solution of 50 mM Tris/1 mM MgCl<sub>2</sub>/50% glycerol, pH 8.0, containing NaN<sub>3</sub> as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL in a stock solution of 50% glycerol/50% PBS, pH 7.4. *No preservative added.* Store at 2-8°C or long-term at -20°C.
- The biotin (BIOT) conjugate is supplied as 1.0 mg in 2.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

## Warning

---

Some reagents contain sodium azide. Please refer to product specific SDS.

## References

---

1. Riedel M, Müller B, Wintermantel E. Protein adsorption and monocyte activation on germanium nanopyrramids. *Biomaterials*. 2001;22:2307-16. (ELISA, Purification)
2. Liu J, Yu J, Wang M, Liu Q, Zhang W, Deng C, et al. Serodiagnosis of Neospora caninum infection in cattle using a recombinant tNcSRS2 protein-based ELISA. *Vet Parasitol*. 2007;143:358-63. (ELISA, WB)
3. Juleff N, Windsor M, Lefevre EA, Gubbins S, Hamblin P, Reid E, et al. Foot-and-mouth disease virus can induce a specific and rapid CD4<sup>+</sup> T-cell-independent neutralizing and isotype class-switched antibody response in naïve cattle. *J Virol*. 2009;83:3626-36. (ELISA)
4. Cui J, Zhang B, Lin Y, Wang Q, Qian A, Nakane A, et al. Protective effect of glutathione S-transferase-fused mutant staphylococcal enterotoxin C against Staphylococcus aureus-induced bovine mastitis. *Vet Immunol Immunopathol*. 2010;135:64-70. (ELISA)
5. Otto W, Najnigier B, Stelmasiak T, Robins-Browne RM. Randomized control trials using a tablet formulation of hyperimmune bovine colostrum to prevent diarrhea caused by enterotoxigenic Escherichia coli in volunteers. *Scand J Gastroenterol*. 2011;46:862-8. (ELISA)
6. Mertens M, Vatansever Z, Mrenoshki S, Krstevski K, Stefanovska J, Djadjovski I, et al. Circulation of Crimean-Congo Hemorrhagic Fever Virus in the former Yugoslav Republic of Macedonia revealed by screening of cattle sera using a novel enzyme-linked immunosorbent assay. *PLoS Negl Trop Dis*. 2015;9(3):e0003519. (ELISA, ICC)
7. Guirnalda P, Murphy NB, Nolan D, Black SJ. Anti-Trypanosoma brucei activity in Cape buffalo serum during the cryptic phase of parasitemia is mediated by antibodies. *Int J Pathol*. 2007;37:1391-9. (FC)
8. Basso W, Lesser M, Grimm F, Hilbe M, Sydler T, Trösch L, et al. Bovine besnoitiosis in Switzerland: imported cases and local transmission. *Vet Parasitol*. 2013;198:265-73. (ICC)
9. Foucras G, Corbière F, Tasca C, Pichereaux C, Caubet C, Trumel C, et al. Alloantibodies against MHC class I: a novel mechanism of neonatal pancytopenia linked to vaccination. *J Immunol*. 2011;187:6564-70. (IP)
10. Odunuga OO, Shazhko A. Ammonium sulfate precipitation combined with liquid chromatography is sufficient for purification of bovine serum albumin that is suitable for most routine laboratory applications. *Bio Chem Comp*. 2013; 1:3 <http://dx.doi.org/10.7243/2052-9341-1-3>. (WB)
11. Indyk HE, Williams JW, Patel HA. Analysis of denaturation of bovine IgG by heat and high pressure using an optical biosensor. *Int Dairy J*. 2008;18:359-66. (Surface Plasmon Resonance)