




## Goat Anti-Mouse Lambda

Cat. No.	Format	Size
1060-01	Purified (UNLB)	1.0 mg
1060-02	Fluorescein (FITC)	1.0 mg
1060-03	Rhodamine (TRITC)	1.0 mg
1060-04	Alkaline Phosphatase (AP)	1.0 mL
1060-05	Horseradish Peroxidase (HRP)	1.0 mL
1060-07	Texas Red® (TXRD)	1.0 mg
1060-08	Biotin (BIOT)	1.0 mg
1060-09	R-phycoerythrin (PE)	0.5 mg
1060-32	Alexa Fluor® 555 (AF555)	1.0 mg

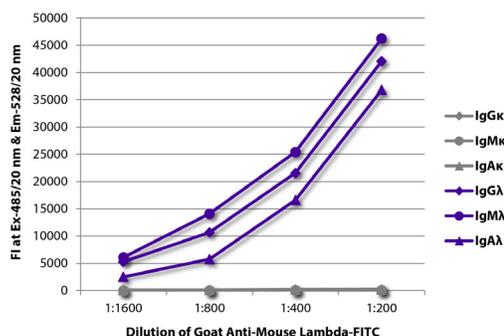

FLISA plate was coated with purified mouse IgGκ, IgMκ, IgAκ, IgGλ, IgMλ, and IgAλ. Immunoglobulins were detected with serially diluted Goat Anti-Mouse Lambda-FITC (SB Cat. No. 1060-02).

### Description

<b>Specificity</b>	Reacts with mouse λ light chains
<b>Source</b>	Pooled antisera from goats hyperimmunized with mouse λ light chains
<b>Cross Adsorption</b>	Mouse κ light chains; may react with λ light chains from other species
<b>Purification</b>	Affinity chromatography on mouse λ light chains covalently linked to agarose

### Applications

Quality tested applications include –

ELISA<sup>1-6,23</sup>  
 FLISA  
 FC<sup>3,5,7-13,24</sup>

Other referenced applications include –

ELISPOT<sup>6,7</sup>  
 IHC-FS<sup>2,13-16</sup>  
 IHC-PS<sup>17,23</sup>  
 ICC<sup>1,4,18</sup>  
 WB<sup>19,20,25</sup>  
 IP<sup>21,22</sup>  
 Stim<sup>23,24</sup>  
 Sep<sup>25</sup>

### Working Dilutions

<b>ELISA</b>	AP conjugate	1:2,000 – 1:4,000
	HRP conjugate	1:4,000 – 1:8,000
	BIOT conjugate	1:5,000 – 1:20,000
<b>FLISA</b>	FITC, TRITC, TXRD, and AF555 conjugates	1:100 – 1:400
	PE conjugate	≤ 1 µg/mL
<b>Flow Cytometry</b>	FITC and BIOT conjugates	≤ 1 µg/10 <sup>6</sup> cells
	PE conjugate	≤ 0.1 µg/10 <sup>6</sup> cells
	For flow cytometry, the suggested use of these reagents is in a final volume of 100 µL	
<b>Other Applications</b>	Since applications vary, you should determine the optimum working dilution for the product that is appropriate for your specific need.	

**For Research Use Only. Not for Diagnostic or Therapeutic Use.**

## Handling and Storage

- The purified (UNLB) antibody is supplied as 1.0 mg purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. *No preservatives or amine-containing buffer salts added.* Store at 2-8°C.
- The fluorescein (FITC), rhodamine (TRITC), and Texas Red® (TXRD) conjugates are supplied as 1.0 mg in 1.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL in a stock solution of 50 mM Tris/1 mM MgCl<sub>2</sub>/50% glycerol, pH 8.0, containing NaN<sub>3</sub> as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL in a stock solution of 50% glycerol/50% PBS, pH 7.4. *No preservative added.* Store at 2-8°C or long-term at -20°C.
- The biotin (BIOT) conjugate is supplied as 1.0 mg in 2.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- The R-phycoerythrin (PE) conjugate is supplied as 0.5 mg in 1.0 mL of PBS/NaN<sub>3</sub> and a stabilizing agent. Store at 2-8°C. **Do not freeze!**
- The Alexa Fluor® 555 (AF555) conjugate is supplied as 1.0 mg in 1.0 mL of PBS/NaN<sub>3</sub>. Store at 2-8°C.
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

## Warning

Some reagents contain sodium azide. Please refer to product specific SDS.

## References

1. Gollahon KA, Hagman J, Brinster RL, Storb U. Ig λ-producing B cells do not show feedback inhibition of gene rearrangement. *J Immunol.* 1988;141:2771-80. (ELISA, ICC)
2. Mandik-Nayak L, Nayak S, Sokol C, Eaton-Bassiri A, Madaio MP, Caton AJ, et al. The origin of anti-nuclear antibodies in bcl-2 transgenic mice. *Int Immunol.* 2000;12:353-64. (ELISA, IHC-FS)
3. Shivtiel S, Leider N, Melamed D. Receptor editing in CD45-deficient immature B cells. *Eur J Immunol.* 2002;32:2264-73. (ELISA, FC)
4. Ghumra A, Semblat J, McIntosh RS, Raza A, Rasmussen IB, Braathen R, et al. Identification of residues in the C<sub>μ</sub>4 domain of polymeric IgM essential for interaction with Plasmodium falciparum erythrocyte membrane protein 1 (PfEMP1). *J Immunol.* 2008;181:1988-2000. (ELISA, ICC)
5. Liu X, Shen S, Manser T. Influence of B cell antigen receptor expression level on pathways of B cell tolerance induction. *J Immunol.* 2009;182:398-407. (ELISA, FC)
6. Willems KN. Regulation of humoral immunity by Pim kinases [dissertation]. Worcester (USA): University of Massachusetts; 2011. (ELISA, ELISPOT)
7. Seagal J, Edry E, Keren Z, Leider N, Benny O, Machluf M, et al. A fail-safe mechanism for negative selection of isotype-switched B cell precursors is regulated by the Fas/FasL pathway. *J Exp Med.* 2003;198:1609-19. (ELISPOT, FC)
8. Adachi T, Wakabayashi C, Nakayama T, Yakura H, Tsubata T. CD72 negatively regulates signaling through the antigen receptor of B cells. *J Immunol.* 2000;164:1223-9. (FC)
9. Su Y, Flemming A, Wossning T, Hobeika E, Reth M, Jumaa H. Identification of a pre-BCR lacking surrogate light chain. *J Exp Med.* 2003;198:1699-706. (FC)
10. Hayashi K, Nijima T, Goitsuka R, Kitamura D. Impaired receptor editing in the primary B cell repertoire of BASH-deficient mice. *J Immunol.* 2004;173:5980-8. (FC)
11. Vettermann C, Herrmann K, Albert C, Roth E, Bösl MR, Jäck H. A unique role for the λ5 nonimmunoglobulin tail in early B lymphocyte development. *J Immunol.* 2008;181:3232-42. (FC)
12. Heizmann B, Kastner P, Chan S. Ikaros is absolutely required for pre-B cell differentiation by attenuating IL-7 signals. *J Exp Med.* 2013;210:2823-32. (FC)
13. Aslam M, Kishi Y, Tsubata T. Excess CD40L does not rescue anti-DNA B cells from clonal anergy. *F1000Res.* 2013;2:218. (FC, IHC-FS)
14. Fukuda T, Yoshida T, Okada S, Hatano M, Miki T, Ishibashi K, et al. Disruption of the Bcl6 gene results in an impaired germinal center formation. *J Exp Med.* 1997;186:439-48. (IHC-FS)
15. Mandik-Nayak L, Huang G, Sheehan KC, Erikson J, Chaplin DD. Signaling through TNF receptor p55 in TNF-α-deficient mice alters the CXCL13/CCL19/CCL21 ratio in the spleen and induces maturation and migration of anergic B cells into the B cell follicle. *J Immunol.* 2001;167:1920-8. (IHC-FS, ICC)
16. Rouaud P, Saintamand A, Saad F, Carrion C, Lecardeur S, Cogné M, et al. Elucidation of the enigmatic IgD class-switch recombination via germline deletion of the IgH 3' regulatory region. *J Exp Med.* 2014;211:975-85. (IHC-FS)
17. Weinstein JS, Nacionales DC, Lee PY, Kelly-Scumpia KM, Yan X, Scumpia PO, et al. Colocalization of antigen-specific B and T cells within ectopic lymphoid tissue following immunization with exogenous antigen. *J Immunol.* 2008;181:3259-67. (IHC-PS)
18. Ghumra A, Shi J, McIntosh RS, Rasmussen IB, Braathen R, Johansen F, et al. Structural requirements for the interaction of human IgM and IgA with the human Fcα/μ receptor. *Eur J Immunol.* 2009;39:1147-56. (ICC)
19. Elkabetz Y, Argon Y, Bar-Nun S. Cysteines in C<sub>H</sub>1 underlie retention of unassembled Ig heavy chains. *J Biol Chem.* 2005;280:14402-12. (WB)
20. Sato M, Adachi T, Tsubata T. Augmentation of signaling through BCR containing IgE but not that containing IgA due to lack of CD22-mediated signal regulation. *J Immunol.* 2007;178:2901-7. (WB)
21. Shen Y, Meunier L, Hendershot LM. Identification and characterization of a novel endoplasmic reticulum (ER) DnaJ homologue, which stimulates ATPase activity of BiP in vitro and is induced by ER stress. *J Biol Chem.* 2002;277:15947-56. (IP)
22. Behnke J, Hendershot LM. The large Hsp70 Grp170 binds to unfolded protein substrates in vivo with a regulation distinct from conventional Hsp70s. *J Biol Chem.* 2014;289:2899-907. (IP)
23. Geraldes P, Rebrovich M, Herrmann K, Wong J, Jäck H, Wabl M, et al. Ig heavy chain promotes mature B cell survival in the absence of light chain. *J Immunol.* 2007;179:1659-68. (ELISA, IHC-FS, Stim)
24. Übelhart R, Bach MP, Eschbach C, Wossning T, Reth M, Jumaa H. N-linked glycosylation selectively regulates autonomous precursor BCR function. *Nat Immunol.* 2010;11:759-65. (FC, Stim)
25. Yurovsky VV, Kelsoe G. Pairing of V<sub>H</sub> gene families with the λ1 light chain: evidence for a non-stochastic association. *Eur J Immunol.* 1993;23:1975-9. (WB, Sep)

Texas Red® is a registered trademark of Molecular Probes, Inc.

Alexa Fluor® 488, 647, and 555 are provided under an Intellectual property license from Life Technologies Corporation. The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or otherwise use this product or its components or materials made using this product or its components for Commercial Purposes. Commercial Purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. For information on purchasing a license to this product for any other use, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or [outlicensing@lifetech.com](mailto:outlicensing@lifetech.com).