




Bovine Type II Collagen

Cat. No.	Format	Size
1220-02	Purified Protein - Lyophilized	0.25 mg
1220-02S	Purified Protein - Solution	0.5 ma

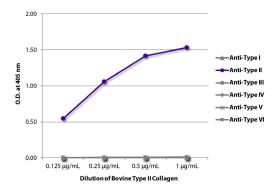

ELISA plate was coated with serially diluted Bovine Type II Collagen (SB Cat. No. 1220-02). Purified collagen was detected with Goat Anti-Type I Collagen-BIOT (SB Cat. No. 1310-08), Goat Anti-Type II Collagen-BIOT (SB Cat. No. 1320-08), Goat Anti-Type III Collagen-BIOT (SB Cat. No. 1330-08), Goat Anti-Type IV Collagen-BIOT (SB Cat. No. 1340-08), Goat Anti-Type V Collagen-BIOT (SB Cat. No. 1350-08), and Goat Anti-Type VI Collagen-BIOT (SB Cat. No. 1360-08) followed by Streptavidin-HRP (SB Cat. No. 7100-05).

Overview

Source Cartilage

Purification Controlled and limited pepsin digestion followed by selective salt precipitation

Purity > 90% by SDS-PAGE

Alternate Name(s) COL2A1

Description

Collagen is the main structural protein in the extracellular space and is the most abundant protein in the ECM. Collagens are divided into two classes - fibril (types I, II, III, V) and non-fibril (types IV, VI). Type II collagen is the main component of cartilage and mutations are associated several types of chondrodysplasia. Type II collagen is formed by homotrimers of $\alpha 1(II)$ chains.

Applications

ELISA – Quality tested ^{1,10} SDS-PAGE – Quality tested ² WB – Reported in literature ^{3,4} *In vivo* Assays – Reported in literature ⁵⁻⁸ Coating Material for – ECM Interaction Studies – Reported in literature ⁹ Migration Studies – Reported in literature ¹⁰

Handling and Storage

- The purified protein is supplied as a solution of 0.5 mg collagen in 1.0 mL of 500 mM acetic acid or 0.25 mg collagen lyophilized from 500 mM acetic acid. Store at 2-8°C.
- Reconstitute lyophilized protein in 500 mM acetic acid.
- Reagents are stable for the period shown on the label if stored as directed.

Warning

Reagents contain acetic acid. Please refer to product specific SDS.

For Research Use Only. Not for Diagnostic or Therapeutic Use.

Email: info@southernbiotech.com • Website: www.southernbiotech.com

References

- 1. Frenz DA, Liu W, Williams JD, Hatcher V, Galinovic-Schwartz V, Flanders KC, et al. Induction of chondrogenesis: requirement for synergistic interaction of basic fibroblast growth factor and transforming growth factor-beta. Development. 1994;120:415-24. (ELISA, Standard Curve)
- Geng Y, McQuillan D, Roughley PJ. SLRP interaction can protect collagen fibrils from cleavage by collagenases. Matrix Biol. 2006;25:484-91. (SDS-PAGE)
- 3. Patti JM, Jonsson H, Guss B, Switalski LM, Wiberg K, Lindberg M, et al. Molecular characterization and expression of a gene encoding a Staphylococcus aureus collagen adhesin. J Biol Chem. 1992;267:4766-72. (WB)
- Lee ER, Lamplugh L, Kluczyk B, Leblond CP, Mort JS. Neoepitopes reveal the features of type II collagen cleavage and the identity of a collagenase involved in the transformation of the epiphyses anlagen in development. Dev Dyn. 2009;238:1547-63. (WB)
- 5. Little RG II, Gazzano-Santoro H, Parent JB, inventors; Xoma Corporation, assignee. Therapeutic uses of bactericidal/permeability increasing protein products. United States patent US 1994/5348942A. 1994 Sep 20. (*In vivo* Assays)
- Burlingham WJ, Love RB, Jankowska-Gan E, Haynes LD, Xu Q, Bobadilla JL, et al. IL-17-dependent cellular immunity to collagen type V predisposes to obliterative bronchiolitis in human lung transplants. J Clin Invest. 2007;117:3498-506. (In vivo Assays)
- 7. Iwai K, Ishii M, Ohshima S, Miyatake K, Saeki Y. Abundant expression of tetraspanin CD9 in activated osteoclasts in ovariectomy-induced osteoporosis and in bone erosions of collagen-induced arthritis. Rheumatol Int. 2008;28:225-31. (*In vivo* Assays)
- 8. Dart ML, Jankowska-Gan E, Huang G, Roenneburg DA, Keller MR, Torrealba JR, et al. Interleukin-17-dependent autoimmunity to collagen type V in atherosclerosis. Circ Res. 2010;107:1106-16. (*In vivo* Assays)
- 9. Hocking AM, Strugnell RA, Ramamurthy P, McQuillan DJ. Éukaryotic expression of recombinant biglycan. Post-translational processing and the importance of secondary structure for biological activity. J Biol Chem. 1996;271:19571-7. (Coating, ECM Interaction Studies)
- Hidaka C, Cheng C, Alexandre D, Bhargava M, Torzilli PA. Maturational differences in superficial and deep zone articular chondrocytes. Cell Tissue Res. 2006;323:127-35. (ELISA, Standard Curve, Coating, Migration Studies)

Tel: 205.945.1774 • U.S. and Canada: 800.722.2255 • Fax: 205.945.8768 Email: info@southernbiotech.com • Website: www.southernbiotech.com