




Mouse Anti-Human IgG₂ Fd

Cat. No.	Format	Size
9080-01	Purified (UNLB)	0.5 mg
9080-02	Fluorescein (FITC)	0.5 mg
9080-04	Alkaline Phosphatase (AP)	1.0 mL
9080-05	Horseradish Peroxidase (HRP)	1.0 mL
9080-08	Biotin (BIOT)	0.5 mg
9080-09	R-phycoerythrin (PE)	0.1 mg
9080-30	Alexa Fluor® 488 (AF488)	0.1 mg
9080-31	Alexa Fluor® 647 (AF647)	0.1 mg
9080-32	Alexa Fluor® 555 (AF555)	0.1 mg

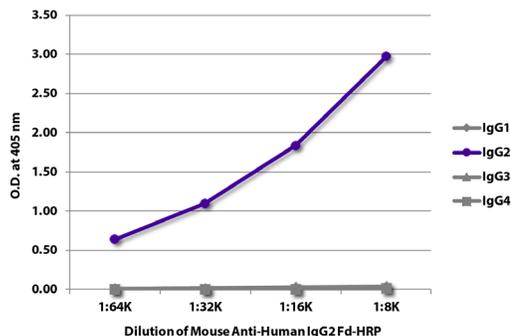

ELISA plate was coated with purified human IgG₁, IgG₂, IgG₃, and IgG₄. Immunoglobulins were detected with serially diluted Mouse Anti-Human IgG₂ Fd-HRP (SB Cat. No. 9080-05).

Overview

Clone	HP6014
Isotype	Mouse (BALB/c) IgG _{1κ}
Immunogen	Human IgG ₂ myeloma protein
Specificity	Human IgG ₂ Fd; Mr 146 kDa

Applications

ELISA – Quality tested ²⁻¹²
 FLISA – Quality tested
 FC – Reported in literature ¹⁹
 IHC-FS – Reported in literature ¹³⁻¹⁵
 IHC-PS – Reported in literature ^{12,16,17}
 WB – Reported in literature ^{13,18,19}
 IP – Reported in literature ¹
 Purification – Reported in literature ¹⁹
 Depletion – Reported in literature ²⁰

Working Dilutions

ELISA	AP conjugate	1:500 – 1:1,000
	HRP conjugate	1:1,000 – 1:2,000
	BIOT conjugate	1:5,000 – 1:10,000

FLISA	FITC, AF488, and AF555 conjugates	1:100 – 1:400
	PE and AF647 conjugates	≤ 1 µg/mL

Other Applications Since applications vary, you should determine the optimum working dilution for the product that is appropriate for your specific need.

For Research Use Only. Not for Diagnostic or Therapeutic Use.

Handling and Storage

- The purified (UNLB) antibody is supplied as 0.5 mg of purified immunoglobulin in 1.0 mL of borate buffered saline, pH 8.2. *No preservatives or amine-containing buffer salts added.* Store at 2-8°C.
- The fluorescein (FITC) conjugate is supplied as 0.5 mg in 1.0 mL of PBS/NaN₃. Store at 2-8°C.
- The alkaline phosphatase (AP) conjugate is supplied as 1.0 mL in a stock solution of 50 mM Tris/1 mM MgCl₂/50% glycerol, pH 8.0, containing NaN₃ as preservative. Store at 2-8°C or long-term at -20°C.
- The horseradish peroxidase (HRP) conjugate is supplied as 1.0 mL in a stock solution of 50% glycerol/50% PBS, pH 7.4. No preservative added. Store at 2-8°C or long-term at -20°C.
- The biotin (BIOT) conjugate is supplied as 0.5 mg in 1.0 mL of PBS/NaN₃. Store at 2-8°C.
- The R-phycoerythrin (PE) conjugate is supplied as 0.1 mg in 1.0 mL of PBS/NaN₃ and a stabilizing agent. Store at 2-8°C. **Do not freeze!**
- The Alexa Fluor® 488 (AF488), Alexa Fluor® 555 (AF555), and Alexa Fluor® 647 (AF647) conjugates are supplied as 0.1 mg in 0.2 mL of PBS/NaN₃. Store at 2-8°C.
- Protect fluorochrome-conjugated forms from light. Reagents are stable for the period shown on the label if stored as directed.

Warning

Some reagents contain sodium azide. Please refer to product specific SDS.

References

1. Reimer CB, Phillips DJ, Aloisio CH, Moore DD, Galland GG, Wells TW, et al. Evaluation of thirty-one mouse monoclonal antibodies to human IgG epitopes. *Hybridoma*. 1984;3:263-75. (Immunogen, IP)
2. Cui Z, Wang H, Zhao M. Natural autoantibodies against glomerular basement membrane exist in normal human sera. *Kidney Int*. 2006;69:894-9. (ELISA)
3. Gao Y, Ye H, Yu F, Guo X, Zhao M. Anti-myeloperoxidase IgG subclass distribution and avidity in sera from patients with propylthiouracil-induced antineutrophil cytoplasmic antibodies associated vasculitis. *Clin Immunol*. 2005;117:87-93. (ELISA)
4. Chapuy-Regaud S, Nogueira L, Clavel C, Sebbag M, Vincent C, Serre G. IgG subclass distribution of the rheumatoid arthritis-specific autoantibodies to citrullinated fibrin. *Clin Exp Immunol*. 2005;139:542-50. (ELISA)
5. Liu L, Chen M, Yu F, Zhao M, Wang H. IgG subclass distribution, affinity of anti-myeloperoxidase antibodies in sera from patients with Wegener's granulomatosis and microscopic polyangiitis. *Nephrology*. 2008;13:629-35. (ELISA)
6. Xie L, Gao Y, Li M, Lu G, Guo X. Distribution of immunoglobulin G subclasses of anti-thyroid peroxidase antibody in sera from patients with Hashimoto's thyroiditis with different thyroid functional status. *Clin Exp Immunol*. 2008;154:172-6. (ELISA)
7. Zhao J, Yan Y, Cui Z, Yang R, Zhao M. The immunoglobulin G subclass distribution of anti-GBM autoantibodies against rH α 3(IV)NC1 is associated with disease severity. *Hum Immunol*. 2009;70:425-9. (ELISA)
8. Xu P, Chen M, Cui Z, Zhao M. Influence of myeloperoxidase by anti-myeloperoxidase antibodies and its association with the disease activity in microscopic polyangiitis. *Rheumatology*. 2010;49:2068-75. (ELISA)
9. Xu P, Cui Z, Chen M, Hellmark T, Zhao M. Comparison of characteristics of natural autoantibodies against myeloperoxidase and anti-myeloperoxidase autoantibodies from patients with microscopic polyangiitis. *Rheumatology*. 2011;50:1236-43. (ELISA)
10. Xu P, Chen M, Zhao M. High potential to reverse the inhibition of myeloperoxidase by ceruloplasmin of anti-myeloperoxidase autoantibodies of IgG3 subclass. *Autoimmunity*. 2012;45:218-25. (ELISA)
11. Shao C, Huo N, Zhao L, Gao Y, Fan X, Zheng Y, et al. The presence of thyroid peroxidase antibody of IgG2 subclass is a risk factor for thyroid dysfunction in chronic hepatitis C patients. *Eur J Endocrinol*. 2013;168:717-22. (ELISA)
12. Jia M, Hu S, Chen J, Qu Z, Liu G, Cui Z, et al. The clinical and immunological features of patients with combined anti-glomerular basement membrane disease and membranous nephropathy. *Kidney Int*. 2014;85:945-52. (ELISA, IHC-PS)
13. Mongan LC, Ockleford CD. Behaviour of two IgG subclasses in transport of immunoglobulin across the human placenta. *J Anat*. 1996;188:43-51. (IHC-FS, WB)
14. Lin J, Markowitz GS, Valeri AM, Kambham N, Sherman WH, Appel GB, et al. Renal monoclonal immunoglobulin deposition disease: the disease spectrum. *J Am Soc Nephrol*. 2001;12:1482-92. (IHC-FS)
15. Greenlee JE, Boyden JW, Pingree M, Brashear HR, Clawson SA, Keeney PM. Antibody types and IgG subclasses in paraneoplastic neurological syndromes. *J Neurol Sci*. 2001;184:131-7. (IHC-FS)
16. Qu Z, Liu G, Li J, Wu L, Tan Y, Zheng X, et al. Absence of glomerular IgG4 deposition in patients with membranous nephropathy may indicate malignancy. *Nephrol Dial Transplant*. 2012;27:1931-7. (IHC-PS)
17. Qu Z, Cui Z, Liu G, Zhao M. The distribution of IgG subclass deposition on renal tissues from patients with anti-glomerular basement membrane disease. *BMC Immunol*. 2013;14:19. (IHC-PS)
18. Harada S, Hata S, Kosada Y, Kondo E. Identification of epitopes recognized by a panel of six anti-human IgG2 monoclonal antibodies. *J Immunol Methods*. 1991;141:89-96. (WB)
19. Keen J, Serghides L, Ayi K, Patel SN, Ayisi J, van Eijk A, et al. HIV impairs opsonic phagocytic clearance of pregnancy-associated malaria parasites. *PLoS Med*. 2007;4(5):e181. (WB, FC, Purification)
20. Forthal DN, Landucci G, Ding H, Kappes JC, Wang A, Thung I, et al. IgG2 inhibits HIV-1 internalization by monocytes, and IgG subclass binding is affected by gp120 glycosylation. *AIDS*. 2011;25:2099-104. (Depletion)

Alexa Fluor® 488, 647, and 555 are provided under an Intellectual property license from Life Technologies Corporation. The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or otherwise use this product or its components or materials made using this product or its components for Commercial Purposes. Commercial Purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. For information on purchasing a license to this product for any other use, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com.

TB9080
05-Dec-21

Corporate Offices: 160 Oxmoor Blvd • Birmingham, AL 35209 • USA **Mailing Address:** P.O. Box 26221 • Birmingham, AL 35260 • USA

Tel: 205.945.1774 • U.S. and Canada: 800.722.2255 • **Fax:** 205.945.8768

Email: info@southernbiotech.com • **Website:** www.southernbiotech.com